




Orbi-Blotter™ Low Speed Orbital Shaker

ORBi-BLOTTER™

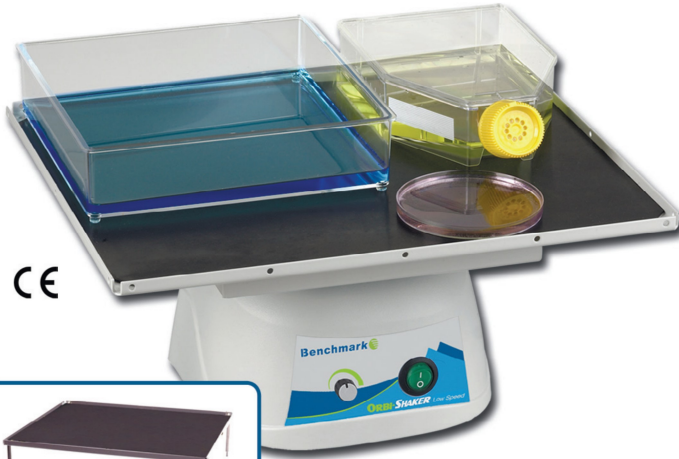

CE

- + Gentle orbital motion for gel applications
- + Large 14 x 12 in. work surface
- + Incubator and cold room safe
- + Variable speed control, 3 to 70 rpm



Optional Stacking Platform (Clearance 3")

The Orbi-Blotter™ features a gentle orbital motion with speed control down to as little as 3 rpm, making it the ideal orbital choice for low speed applications, such as blotting, washing, staining/destaining, etc. A large 14x12 in. work surface is included along with a non-slip rubber mat for holding trays, dishes, plates and other flat vessels.

An optional stacking platform can be purchased as an accessory, effectively doubling the usable work surface. The maintenance free, brushless motor is capable of loads up to 4kg and is safe for use in cold rooms and incubators.

Technical Data:

Speed Range:	Variable, 3 to 70 rpm
Motion / Orbit:	Orbital / 19mm
Load Capacity:	3 x 2L or 9 lbs (4kg)
Platform Size:	14 x 12 in. / 35 x 30 cm
Operating Temp. Range:	+4°C to +45°C
Dimensions: (WxDxH)	14 x 12 x 8 in. 35.5 x 30.5 x 20 cm
Stacking Clearance:	3 in. / 7.6 cm
Weight:	7.5 lbs (3.4kg)
Electrical:	115 or 230V, 50-60Hz, 50W
Warranty:	2 Years

Ordering Information:

Item No.	Description:
BT30*	Orbi-Blotter™ Low Speed Orbital Shaker
BR2000-STACK	Stacking platform with 14 x 12" flat mat (3" clearance)
BR2000-SP	STACK-PLUS Extender set, adds 1.25" clearance to stacking platform

*115V with US plug. For 230v (EU plug), please add (-E).

Benchmotion™ Shaker Selection App

BenchMotion Shaker Selection App



Available on:



Scan me

Simply download the app on your phone or mobile device, or visit www.benchmotion.com to view a video demonstration showing all of the shaking motions available from Benchmark Scientific. Once the preferred motion is identified, click or touch on the "see details" option to quickly navigate your way to the most appropriate model(s) for your application.