




BenchMixer™ Vortex Mixers

BenchMixer™ V2

- Reversible top: Use “Cup” or “Flat/Dimpled” Surface
- Powerful motor for INSTANT vortexing
- Q-Drive™ dynamic balancing system
- Variable speed with touch or continuous operation

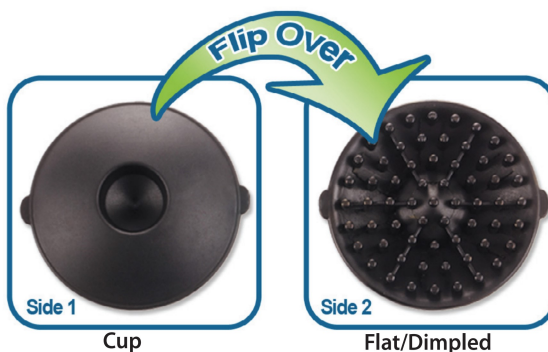


Reversible Top

The BenchMixer™ V2 utilizes a smooth and powerful mixing technology similar to the original BenchMixer, with the advantage of a uniquely designed reversible rubber head that can be attached to the mixer in two orientations. For the vortexing of single tubes, the cupped surface would generally be preferred. However, for applications requiring a flat/dimpled surface for mixing uniquely shaped vessels (or a handful of multiple test-tubes), the head can easily be flipped over and reattached in reverse orientation.

Like the original BenchMixer, the drive system is balanced with Q-Drive™ technology, providing a smooth vortexing action, while minimizing noise and excessive vibration. This results in quieter, more efficient and longer lasting operation. With a maximum speed up to 3,200 rpm, the mixer provides an instant vortex of even the largest sample sizes, including full 50ml tubes.

Available Accessories:

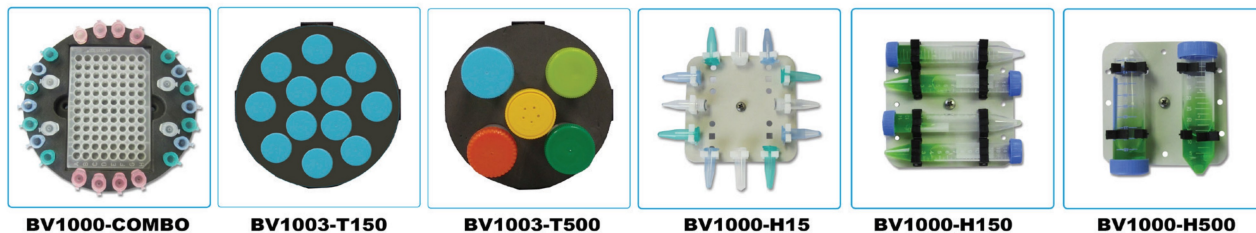

BV1000-COMBO

BV1003-T150

BV1003-T500

BV1000-H15

BV1000-H150

BV1000-H500

Technical Data:

Speed Range:	200 to 3,200 rpm
Operating Modes:	Touch or continuous
Orbit:	3mm
Operating Temp. Range:	+4° C to +45°C
Dimensions: (W x D x H)	5 x 6.3 x 6.7 in.* 13 x 16 x 17 cm
Weight:	8.4 lbs (3.8kg)
Electrical:	115V, 60Hz or 230V, 50Hz/60Hz, 180W
Warranty:	2 Years

*Dimensions of BV1000 model: 5.2 x 6.7 x 6.7 in./13 x 17 x 17 cm

Ordering Information:

Item No.	Description:
BV1003*	BenchMixer V2 Vortexer (with reversible head)
BV1000*	BenchMixer Vortexer (original) with standard cup head
BV1000-FLAT	Optional 3 in. flat head attachment (for BV1000 only)
BV1000-COMBO	Optional combination head for tubes and microplates (38x1.5/2.0ml & 28x0.5ml or one microplate)
BV1003-T150	Optional head for 12 x 15ml tubes (BV1003 only)
BV1003-T500	Optional head for 5 x 50ml tubes (BV1003 only)
BV1000-H15	Optional head for 12 x 1.5/2.0ml, horizontal
BV1000-H150	Optional head for 4 x 15ml tubes, horizontal
BV1000-H500	Optional head for 2 x 50ml tubes, horizontal



BV1003



BV1000

*115V with US plug. For 230V/50Hz with EU plug, please add (-E). For 230V/60Hz with EU plug, please add (-E2).