




- + The only vortexer available with **Multi-Head™**
- + Powerful motor for **INSTANT** vortexing
- + **Q-Drive™** dynamic balancing system
- + Variable speed from 200 to 3,200 rpm



Unique Multi-Head allows for standard vortexing, **PLUS** holds up to 8 microtubes

Available Accessories:

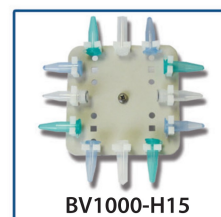

The Mortexer™ is MORE than your standard vortex mixer. The unique design of the Multi-Head™ (included) features a traditional cup head for general purpose mixing, plus holds up to eight microcentrifuge tubes. The drive system is balanced with Q-Drive™ technology, providing a smooth vortexing action, while minimizing noise and excessive vibration. This results in quieter, more efficient and longer lasting operation.

Other features include: “no-walk” suction feet, cold room/incubator compatibility and spill-proof electronic controls.

Technical Data:

Speed Range:	200 to 3,200 rpm
Operating Modes:	Touch or continuous
Orbit:	3mm
Operating Temp. Range:	+4° C to +45° C
Dimensions: (W x D x H)	5 x 6.3 x 6.7 in. 13 x 16 x 17 cm
Weight:	8.4 lbs (3.8kg)
Electrical:	115V, 60Hz or 230V, 50Hz*, 180W
Warranty:	2 Years

Ordering Information:

Item No.	Description:
BV1005*	Mortexer Vortex Mixer (with Multi-Head)
BV1000-FLAT	Optional 3 inch flat head attachment
BV1000-COMBO	Optional combination head for tubes and microplates (38x1.5/2.0ml & 28 x 0.5ml or one microplate)
BV1000-H15	Optional head for 12 x 1.5/2.0ml, horizontal
BV1000-H150	Optional head for 4 x 15ml tubes, horizontal
BV1000-H500	Optional head for 2 x 50ml tubes, horizontal

*115V with US plug. For 230V/50Hz with EU plug, please add (-E). For 230V/60Hz with EU plug, please add (-E2).